



**Presto Stantest Pvt. Ltd.**

I-42, DLF Industrial Area Phase-1, Delhi Mathura Road, Faridabad 121003, Haryana, India  
P : 9210 903 903, +91 129 4272727, 93111 24302 E : info@prestogroup.com

• Faridabad • Sonapat • Kolkata • Mumbai • Pune • Ahmedabad  
• Chennai • Bangalore • Hyderabad

[www.prestogroup.com](http://www.prestogroup.com)



## Specific Gravity Balance Digital




### Model No. SP 202 (G)

**Specific Gravity** of any elastomer is of considerable importance for its identification. Besides it can also be used for providing a rough assessment of its physical properties. It can also be used for finding if the sample and the products are prepared from the same compound or not. Displays direct specific Gravity digitally, No calculations required. Based on the Archimedes principle.

Designed for use in laboratory applications, it is easy to use two-button operation, a high-contrast LCD display.

\* It is must to sink the sample into the DM water to avoid the floating of sample.

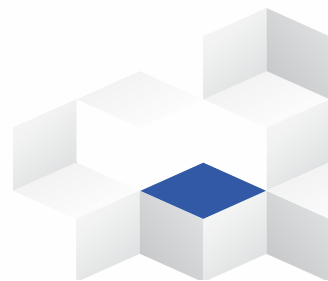



# Presto Stantest Pvt. Ltd.

I-42, DLF Industrial Area Phase-1, Delhi Mathura Road, Faridabad 121003, Haryana, India  
P : 9210 903 903, +91 129 4272727, 93111 24302 E : info@prestogroup.com

- Faridabad • Sonipat • Kolkata • Mumbai • Pune • Ahmedabad
- Chennai • Bangalore • Hyderabad

[www.prestogroup.com](http://www.prestogroup.com)



## Key Specifications:

|                        |  |
|------------------------|--|
| Max. Weighing Capacity | 200 gm   |
| Least Count            | 0.01 gm  |
| Pan size               | 120 mm   |
| Display                | High contrast LCD display                            |
| Calibration            | 200 gm brass weight (Included)                       |
| Linearity              | $\pm 0.01$ g   |
| Specific Gravity Range | 1.0 to 3.0   |
| Medium                 | Lab Grade DM Water                                   |
| Power                  | AC adaptor (Included) or 4 AA battery (not included) |

## Highlights:

|                 |                                |
|-----------------|--------------------------------|
| Hanging Hook    | Yes                            |
| Beaker Platform | Yes                            |
| Supply          | Main Supply & battery operated |

## Optional :

- Glass beaker available on request.

## Feature :

- Superior shock resistance & overload protection
- Software lock switch & integrated transportation lock
- Adjustable footrest with level indicator.

Thank you customers for choosing us as your partners in growth !



\*All the above logos are the sole property of their respective owners and are used purely for depiction purposes only.

## OUR OTHER PRODUCTS :-

|                        |                         |                        |                    |                    |                     |
|------------------------|-------------------------|------------------------|--------------------|--------------------|---------------------|
|                        |                         |                        |                    |                    |                     |
| Box Compression Tester | Scuff Resistance Tester | Tensile Tester Digital | Spectrocolorimeter | Cobb Sizing Tester | Digital GSM Balance |

